

Bendrigg Accessible Equipment



Version 1: March 2026

Accessible equipment

At Bendrigg we understand being away from home can be a challenge. When it comes to packing the car or van you may wonder if it's even possible at all due to all the equipment some individuals need during their stay.

We want to help; hence we offer a range of accessibility equipment for you to borrow during your stay free of charge.

This catalogue is designed to guide you through what is available and the suitability of each item in different locations.


We can't offer training or assessment on each piece of equipment; we recommend only borrowing equipment that is suited to the individual's needs and that you know how to operate yourself.

If you have questions about items in the catalogue, please do call us and we will do our best to assist you.

Pre-booking is required. Some items are very popular, and we may not always have the number of items required. Where this happens, we will do our best to offer you an alternative. If we don't have an alternative, we will ask you to bring your own.



About Wet rooms

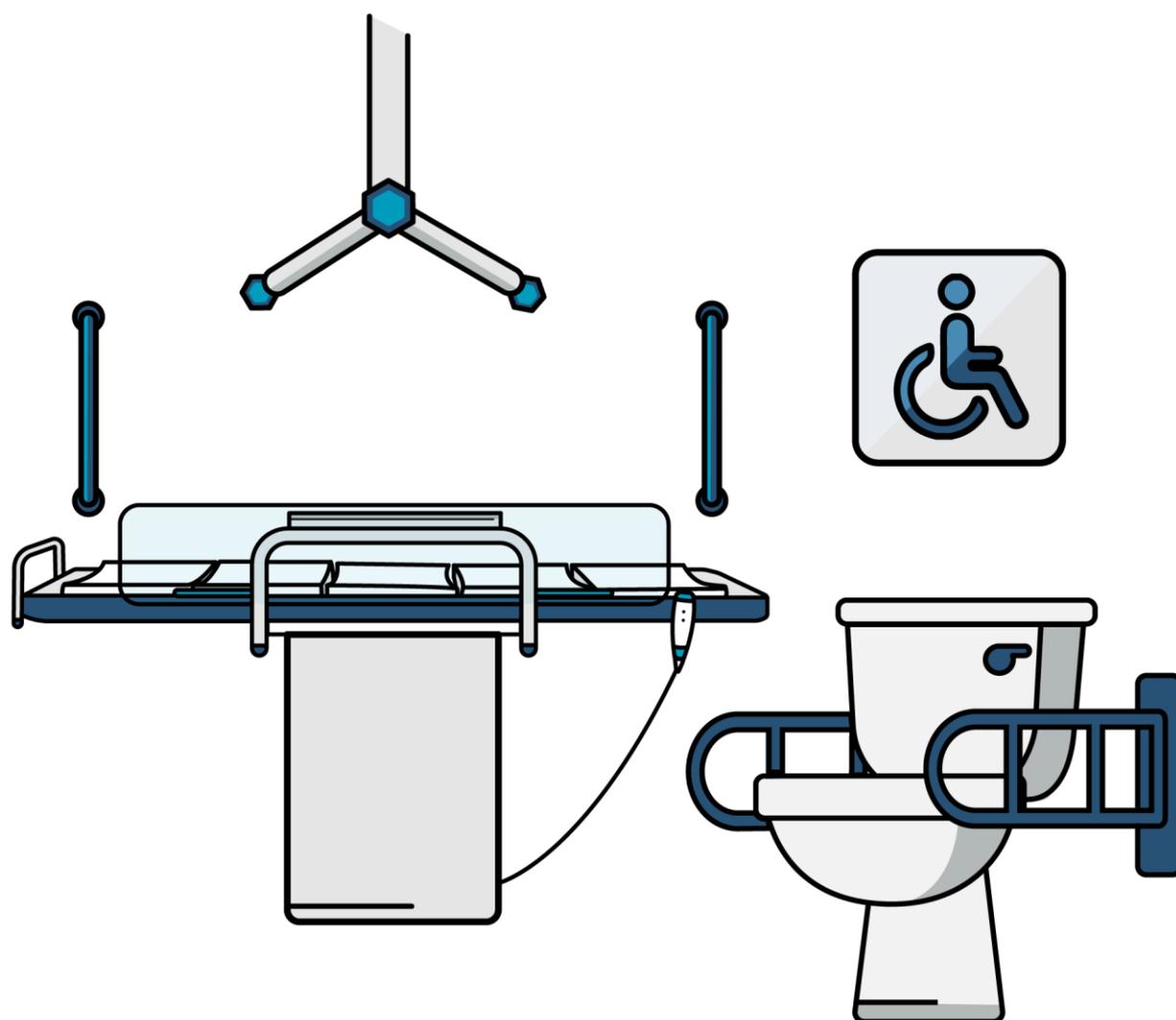

Throughout the Bendrigg site there are a range of wet rooms, all with facilities which we hope will suit your needs.


All hoists throughout the site have loop bars. Please bring your own sling. We have a range of slings for emergency use only; they cannot be booked out.

Large sanitary waste bins are available on request.

Lodge Wet rooms



The lodge is the oldest building; however, it does offer a range of accessibility features within the wet rooms. All wet rooms have automatic lighting, and non-slip flooring.

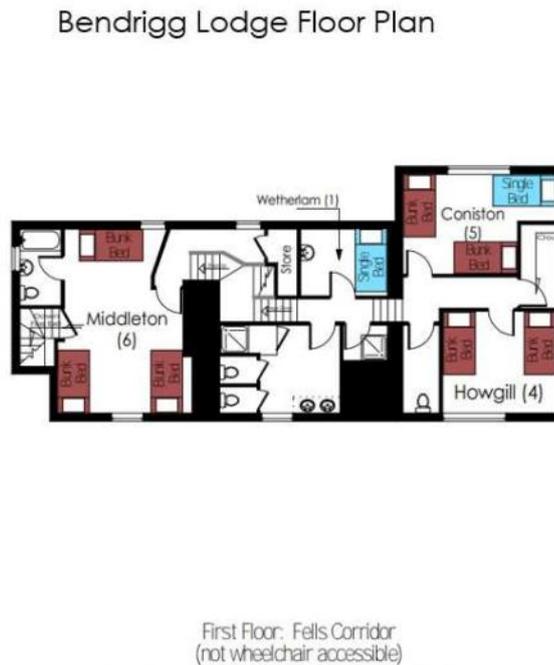

- = Shower Bed
- = Single Track Hoisting System
- = H track hoisting system
- = Profile Bed
- = Single Bed
- = Bunk Bed

NOT TO SCALE: For guidance purposes only

Corridor	Wet room Location	Toilet Height (cm)
Lakes	Opposite Gift Shop	47
Lakes	Opposite Springs Bedroom	Cubicle not accessible for toilet seat.
Lakes	Opposite Lillymere Bedroom	50
Lakes	Opposite Killington Bedroom	47
Clouds	Near Cirrus	45
Clouds	Cumulous En-suite	47
Clouds	Opposite Stratus Bedroom	47
Clouds	Nimbus En-suite	44

Lakes Corridor



On the lakes corridor you will find 4 wet rooms / bathrooms.

The bathroom opposite the gift shop offers:

- Door width - 77cm
- Bath with an adult-sized, wall-mounted drop-down shower bed
- Shelf
- Pedestal sink - 81cm high with wall-mounted grab rails
- Peninsular toilet 47cm high with drop-down grab rails
- Single track hoist over toilet and shower bed. Max. weight 130kg

Opposite Springs bedroom is a wet room / shower room with urinals:

- Door width - 70cm
- 2 shower cubicles with shower curtains (83cm x 100cm)
- 1 small toilet in a cubical, door width 71cm
- 2 urinals 63cm high
- 3 sinks 80cm high
- 2 shelves


Opposite Lillymere bedroom is a wet room:

- Door width – 71cm
- Ceiling H track hoist, max. weight 175kg
- Wall-mounted shower with single grab bar
- Left-sided transfer toilet (50cm height) with drop-down grab bar and wall-mounted grab bars.
- Wall-mounted sink, 75cm high
- Our portable hydraulic shower bed will fit in this wet room, but it is very tight.
- Wall-mounted hand dryer.

Opposite Killington bedroom is a wet room which offers:

- Door width - 80cm
- Wall-mounted shower with single grab bar, and shower curtain
- Right-sided transfer toilet, 47cm high, with drop-down grab bar and wall-mounted grab bars.

Clouds Corridor

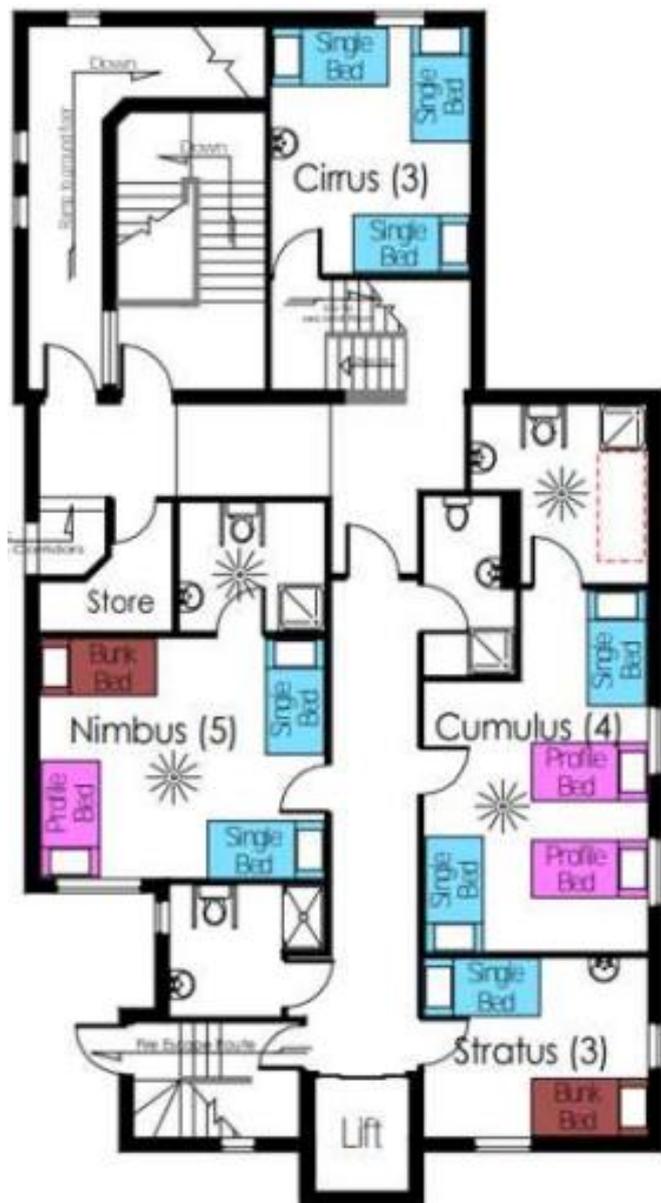

First Floor:
Clouds Corridor
(all accessible except Cirrus)

The Clouds corridor is accessible via ramp, stairs, or the lift.

Wet room immediately through the door on the left (Cirrus' wet room):

- Door width - 88cm
- Left-sided transfer toilet, 45cm high with drop-down grab bar and wall-mounted grab bars.
- Pedestal sink 77cm high with wall-mounted grab bars.
- Wall-mounted height-adjustable shower with grab bars
- Wall-mounted hand dryer.

Cumulus en-suite wet room:

- Door Width - 88cm to Cumulus and 83cm to wet room
- Sink - 78 cm high, with wall-mounted grab bars
- Peninsular toilet 47cm high, with drop-down grab bars
- Wall-mounted shower with grab bars and shower curtain
- Wall-mounted drop-down shower chair with arms, left or front transfer
- Adult-sized, wall-mounted, height-adjustable, profiling, shower bed
- Track hoist travels from bedroom on a single track into the wet room to a junction, which facilitates access to the toilet and shower bed only


Opposite Stratus bedroom the wet room offers:

- Door width - 69cm
- Right-sided transfer toilet 47cm with drop-down grab bar and wall-mounted grab bars
- Pedestal sink - 79cm high, with wall-mounted grab rails
- Shower cubical with curtain (110 x80cm)
- Wall-mounted drop-down shower chair with wall-mounted bars
- Wall-mounted hand dryer

Nimbus en-suite wet room:

- Door width 90cm to Nimbus and 83cm to wet room
- Peninsular toilet height 44cm, with drop-down grab bars
- Wall-mounted shower with grab bars and shower curtain
- Wall-mounted drop-down shower chair with arms, 47cm high
- Track hoist travels from bedroom on a single track into the wet room, which facilitates access to the toilet only
- Our portable hydraulic shower bed will fit in this wet room, but it is very tight

Fells corridor

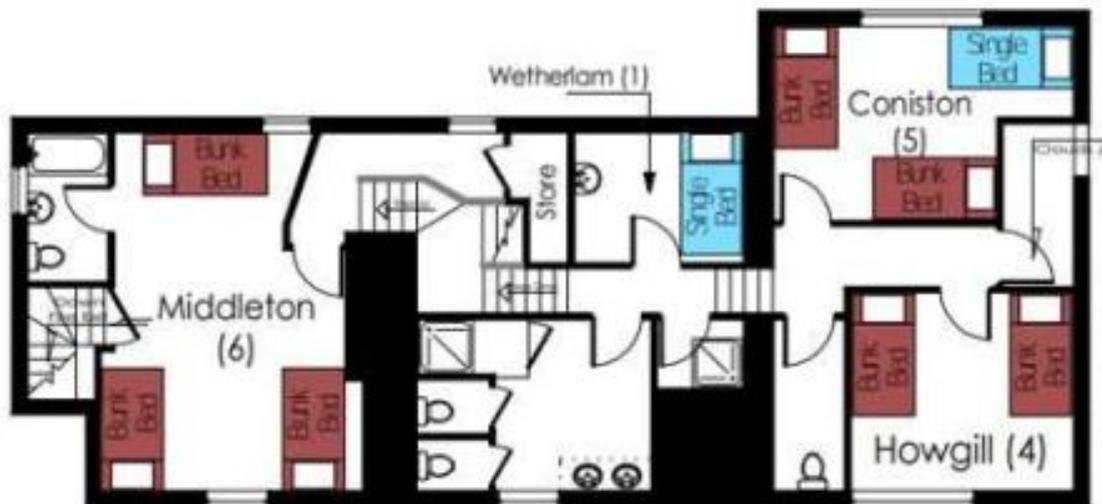

The Fells corridor is accessible via stairs only.

Middleton en-suite bathroom:

- An over-bath shower
- Peninsular toilet
- Pedestal sink

The shower rooms opposite Wetherlam bedroom consist of:

- 2 sinks
- 1 shower cubical
- 2 toilet cubicles
- A separate shower room

Opposite Coniston bedroom:

- Single toilet with sink

Acorn Wet rooms



Acorn House was designed with accessibility and flexibility in mind. All wet rooms have automatic lighting and non-slip flooring.



All Acorn toilets are set at a height of 45cm

The lobby toilet, situated in the entrance hall, is not wheelchair accessible:

- Door width 96cm
- Toilet
- Small sink

Oak and Chestnut wet rooms:

- Door width 96cm to bedroom and 94cm to wet room
- Semi-recessed sink, 85cm high
- Wall-hung toilet, right-sided transfer, 45cm high
- Height-adjustable shower

Hazel, Beech, Silver Birch, and Willow wet rooms:

- Door width 94cm to bedroom and 94cm to wet room
- Wall-mounted sink, 70cm high
- Peninsular toilet, 45cm high, with drop-down grab bars
- Height-adjustable shower + fixed shower
- Wall-mounted drop-down shower chair with arms, 47cm high
- Wall-mounted, fold-down, profiling, powered, height-adjustable shower bed with side; length 155cm width 61cm

Ash wet room:

- Door width 100cm to bedroom and 94cm to wet room
- Wall-mounted sink, 72cm high
- Toilet with drop-down grab bars, right-sided transfer, height 45cm
- Height-adjustable shower + fixed shower
- Wall-mounted drop-down shower chair with arms, 47cm high
- Our portable hydraulic shower bed will fit in this wet room

Cherry and Holly wet rooms:

- Door width 94cm to bedroom and 95cm to wet room
- Semi-recessed sink, 85cm high
- Wall-hung toilet, left-sided transfer, height 45cm
- Height-adjustable shower

Rowan wet room:

- Door width 91cm to bedroom and 93cm to wet room
- Wall-mounted sink, 71cm high
- Toilet with drop-down grab bars, left-sided transfer, height 45cm
- Height-adjustable shower + fixed shower
- Wall-mounted drop-down shower chair with arms, 46cm high
- Our portable hydraulic shower bed will fit in this wet room

Loan equipment


Wet room equipment

Sleeping

Moving About

Comfort and Support

Wet room equipment

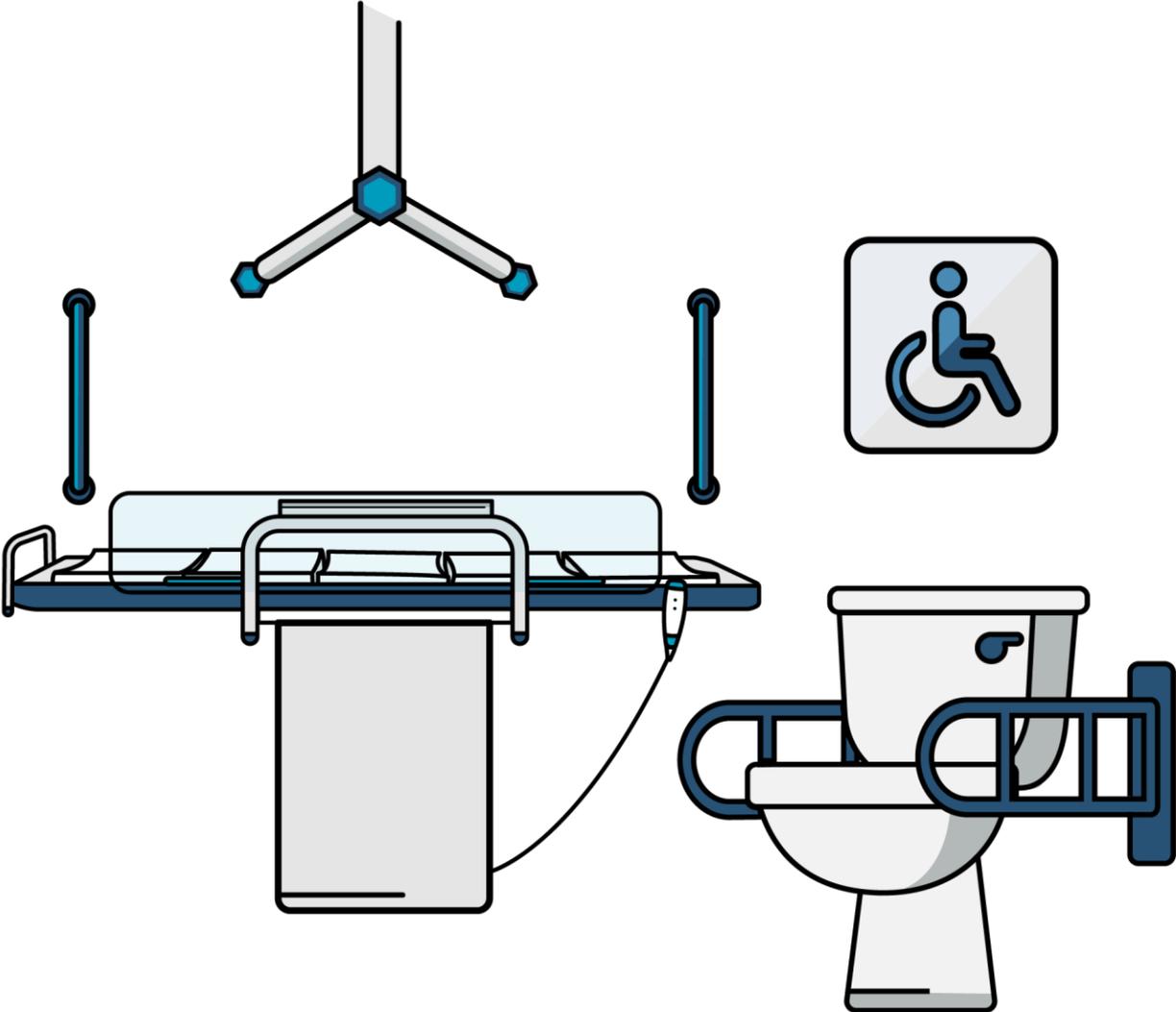

Toilet seat topper (x2)



Summary: raised toilet seat. 100mm raise with side-mounted brackets.

Max Weight Limit: 190kg (30 stone)

Important Info: Large aperture and front cutaway

Cleaning: Wipe clean

Room Use: Not Lilymere or Stratus wet room

Folding commode (x1)



Summary: this folding commode chair provides easy storage and transportation.

Max. user weight limit: 120 kg (18 stone)

Important Info: seat size W38cm, H45cm

Cleaning: wipe down

Room Use: all (with commode bucket). Does not go over toilet.

Freeway (1x 44cm and 1x 55cm)



Summary: the Freeway T40 chair is a highly versatile assistant-propelled shower, toilet, or commode chair.

Max. Weight Limit: 200kg (31 stone)

Special features: padded seat base, back and footrests

Cleaning: wipe down

Room Use: goes over all toilets in Acorn House. Fits over all Clouds and Lakes Corridor toilets except Lilymere.

Freeway self-propelling T70 shower chair (x1)



Summary: can be used as a commode (will not fit over any toilets) or shower chair

Max. Weight Limit: 120kg (18 stone). **Seat:** 45cm wide, 43cm high

Cleaning: wipe clean

Special Features:

- removable armrests
- swing away footrests
- padded seat for extra comfort

Room Use: any wet rooms in Acorn House, and any wet rooms in the Clouds or Lakes corridors (shower chair or commode use only, will not fit over any toilets)

Drive mobile shower commode chair (x1)



Summary: waterproof padded backrest, armrests and seat, all made from soft, easy to clean polyurethane. Seat of the shower commode chair has a cut-out section and commode pan.

Age range: large teen to adult

Max. Weight Limit: 130kg

Important Info: full footplate

Cleaning: wipe clean

Special features: full foot plate

Accessories:

- Commode pan
- Full seat pad.

Room use: Commode use only. Will not fit over any toilets.

Blue shower commode chair (x1)



Summary: Use over a toilet or as a mobile, freestanding toilet chair or in the shower.

Age range: Large teen to adult

Max. Weight Limit: 130kg

Cleaning: wipe down

Accessories:

- Split single footplates
- Soft back support
- Padded arm rests

Room use: Use over all toilets in Acorn. Fits over all Clouds and Lakes Corridor toilets except wet room opposite Lillymere.

Etac toilet swift – Mobil tilt (x1)



Summary: Height-adjustable shower and toilet chair. Soft mesh backrest moulds to fit the user. **We do not have footrests!**

Max. Weight Limit: 160kg (25 stone)

Important Info: detachable and replaceable seat.

Cleaning: wipe clean

Special Features:

- The water-repellent nylon backrest moulds to the user's back for comfort, won't absorb soiled water and dries quickly after use
- Height-adjustable
- Tilt-in-space
- Safe and comfortable seat

Room Use: All except Fells corridor.

Rifton HTS Hygiene Toilet System (2x large, 1x medium)



Summary: Fully adaptable toilet, or can be used as a mobile, freestanding toilet chair or in the shower. Tilt-in-space.

Height-adjustable.

Size adjustable to meet the needs of children and adults with physical disabilities.

Max. user weight: Medium - 68kg. Large - 113kg.

Important Info: headrest, laterals, chest harness, pelvis strap, hip guides, abductor footboard, splash guard and tray.

Cleaning: wipe clean

Room Use: All wet rooms except those on Fells Corridor.

Seahorse (2 x large)



Summary: Use over a toilet or as a mobile, freestanding toilet chair, or in the shower.

Age range: 5 years to adult

Max. Weight Limit: 102kg

Cleaning: Wipe clean

Special features:

- Tilt-in-space
- Height-adjustable head and footrest
- Adjustable 5-point harness
- Removable pommel/splash guard
- Removable seat pan

Room Use: All wet rooms except those on Fells Corridor.

Mobile shower chair (x2)



Summary: This fixed-height, wheeled chair with armrests ensures a secure and comfortable showering experience for users. Designed for use in the shower.

Age range: None specified

Max Weight Limit: 140 kg (22 stone)

Cleaning: Wipe clean

Room use: All wet rooms except those on Fells Corridor.

Folding wall-mounted shower chair (x8)



Summary: Low-profile folding shower seat with legs, back rest and armrests. When not in use, it folds flat against the wall.

Drainage holes in the moulded seat mean there is no slippery standing water.

Max. Weight Limit: 160kg

Cleaning: Wipe down

Size: Seat width 46cm, seat depth 38cm

Room use: Does not need to be requested. Installed in Lodge – Nimbus and Cumulus (Lodge), and Willow, Holly, Silver Birch, Beech, Ash and Rowan (Acorn).

Big Splashy (x1)



Summary: A portable bath seat designed to sit low in the bath to allow children and teenagers to benefit from being immersed in water.

Age range: 7- 14 years

Max. Weight Limit: 55kg

Important Info: Can be used with or without the base.

Cleaning: Wipe clean

Special features:

- Floor sitter has 26 recline positions, from upright (106°) to right back (140°)
- 5 point harness
- Height-adjustable headrest and laterals

Room Use: All wet rooms except those on the Fells corridor (goes on floor or in bath). Can be used as a floor sitter or adaptive seating in other areas.

Wall-mounted shower bed / changing plinth (x3)



Summary: Fixed-height foldable shower bed / changing plinth. When not in use, it folds flat against the wall.

Max. Weight Limit: 20 stone (130kg)

Cleaning: Wipe down

Room use: Does not need to be requested. Installed in Bathroom opposite tuckshop (Lodge), Games room wet room, and wet room at the top of the climbing wall.

Hydraulic shower trolley (x2)



Summary: This mobile, hydraulically operated shower trolley is designed for use in and to provide transport to and from the showering area.

Age range: none specified

Max. Weight Limit: 200kg (31 Stone)


Cleaning: Wipe down

Special features:

- Easy steer mechanism
- Hydraulic height-adjustment allows the carer to move the trolley up and down from either side
- Convenient drop sides and safety catches
- Will comfortably accommodate a large adult

Room use: Fits in Rowan and Ash wet rooms (Acorn). Nimbus or Lillymere wet rooms (Lodge). Note that Nimbus hoist does not go to position of shower bed.

Ropox height-adjustable changing/shower Bed (x5)

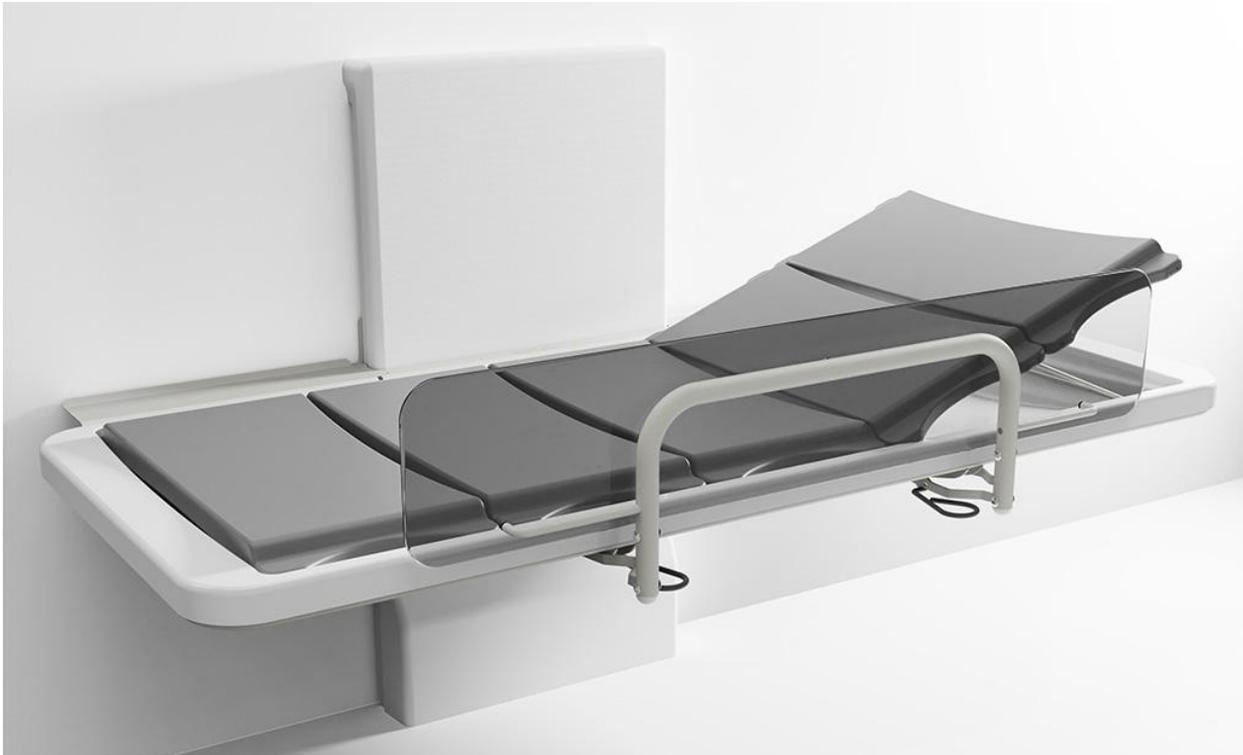

Summary: Wall-mounted and foldable changing and/or shower bed with electric 70 cm height adjustment and back support.

Age range: all

Max. Weight Limit: 220kg

Cleaning: Wipe clean

Room use: Does not need to be requested. Installed in Beech, Hazel, Willow, Silver Birch (Acorn) and Cumulus (Lodge).

Sleeping



BED TIME

Profile beds (x18)



Summary: All powered beds in the Lodge and Acorn height adjust, profile and tilt. They are all fitted with anti-entrapment mesh-padded bed rails so that there is no gap between the mattress platform and side rails.

Important Info: These beds do not need to be requested and cannot be moved into different rooms. Can be used as a standard bed with the rails down.

Room use: Does not need to be requested.

Acorn: Oak (1), Beech (1), Hazel (2), Ash (2), Rowan (2), Willow (2), Silver Birch (2)

Lodge: Springs (2), Lilymere (1), Nimbus (1), Cumulus (2)

Extra-long bed extension & rails for profiling bed (x1)



Summary: For people who are very tall and need to use a profile bed with extra bed length, we have a bed extension pack.

Important Info: Currently fitted to a bed in Willow. Could be moved, if absolutely necessary, to any bed in Acorn.

Special Features: Includes extra long mesh and padded sides.

Cloud cuddle (x2)



Summary: This is a bed tent which is specially developed to allow a child to sleep safely anywhere. Cloud Cuddle is a simple modification that transforms an 'ordinary single' bed into a safe place.

Height Sitting: Max 75cm **Max. User Weight:** 50kg

Cleaning: The zipped-in base sheet can be removed and washed/disinfected separately at 60 degrees Celsius (6 hours drying time) without taking down the Cloud Cuddle.

Special Features: Can be fitted to all divans as well as profile beds as long as sides are not used. Can also be used on a mattress on the floor.

The Safety Sleeper (x1)



Summary: Enclosed bed system designed to keep a child safe inside. Provides a sense of calm and security, and discourages night-wandering. Can be placed on a divan bed or on the floor.

Cleaning: Machine washable & dryable (low heat) (6 hours drying time)

Special Features: One-sided entry and all standard safety features. Built from durable materials and enhanced joints to stand up against heavy use and play

Room Use:

Acorn: If placing on the floor, Oak and Cherry are preferred rooms for use. Also works in Rowan and Ash, but may block interconnecting door and may need to take profile beds into account. Tight fit in Chestnut and Holly.

Lodge: Cumulus, Cirrus and Nimbus. Tight fit in Stratus, Lillymere and Springs.

Safespace® (x1)



Summary: A soft, flexible sleep and relaxation space.

The Safespace® is designed for individuals who require safe surroundings. The flexible sides and soft floor reduce the risk of injury from hard surfaces, providing peace of mind for parents and carers. Our Safespace has a single pedestrian door on one side, and a mesh roof.

Cleaning: Wipe clean

Special Features: Safe, strong and durable for adults and children

Designed to meet complex needs

Significantly reduces the risk of injury from hard surfaces

Calming low stimulation environment

Easily transformed into a sensory space

Room use: Lodge - Springs only. Acorn - Oak and Cherry only, or Games room as a break out space not a bedroom.

Travel cot (x3) (only 1 mattress)



Summary: Travel cot. Suitable from birth.


Max. Weight Limit: up to 15kg

Important Info:

- Mattress dimensions: L97cm, W65cm.

Cleaning: Wipe clean

Moving About

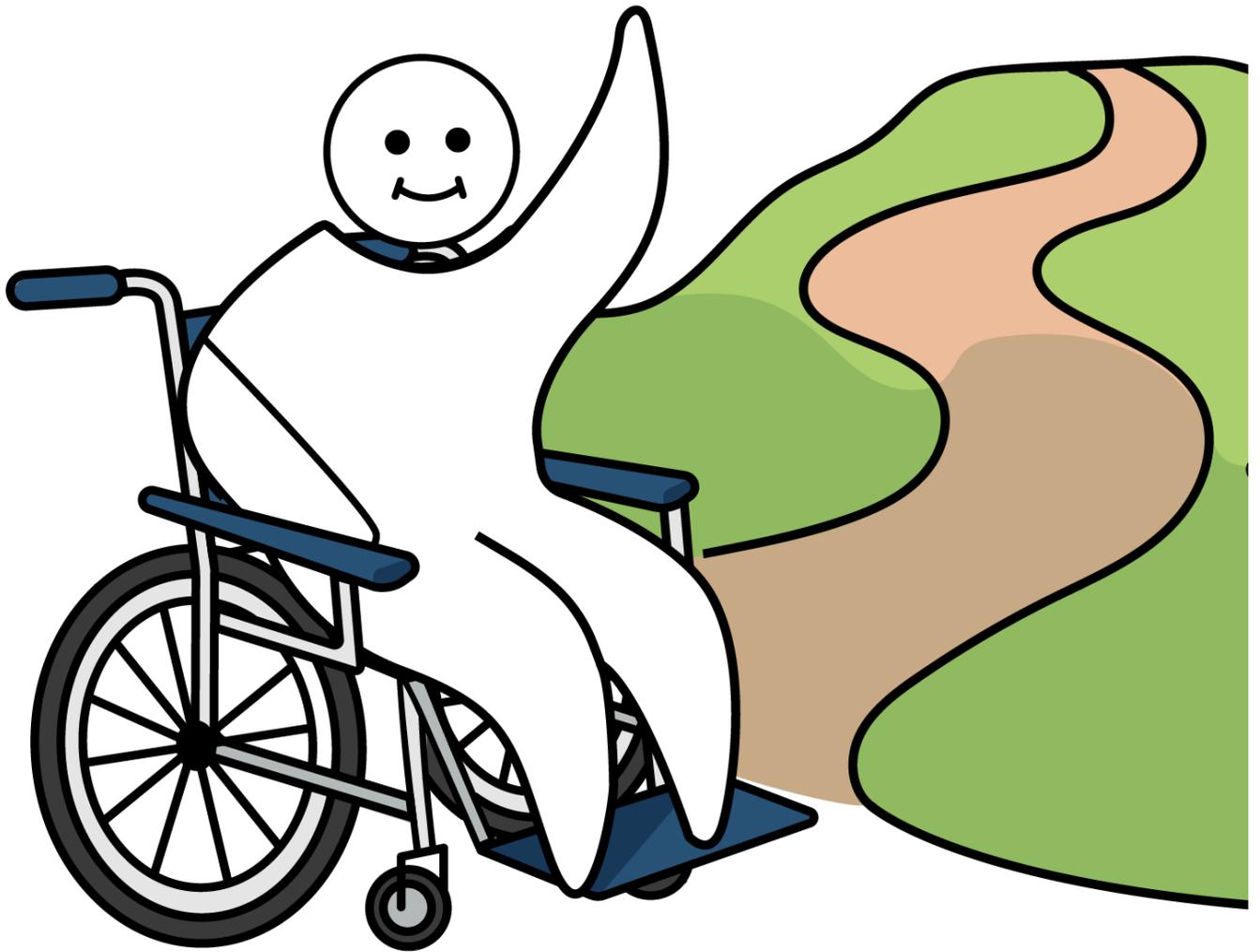

Advance Mobility Pushchairs (x2)



Summary: 3-wheeled pushchair with 5-point harness and height-adjustable foot rest.

Age range: 5-10 years

Max. Weight Limit: 45kg

Cleaning: Washing Machine fabrics 30 deg C. Wipe down frame.

Special features:

Can be used with grey size 2 Special Tomato Soft Touch sitter inserted (shown in picture above).

Accessories:

- Rain cover
- Sun canopy
- Storage basket
- 5-point harness

Sara Stedy (x1)



Summary: The Sara Stedy enables people to take an active role in standing and promote mobility. NOTE: very similar to a QuickMove.

For use by: Those who are 149-193 cm tall (4'11"-6'6")

Max. Weight Limit: 182kg

Cleaning: Wipe clean

Special features:

- The pivoting seat can be moved out to enable the patient to stand and provides an angled seated support during the transfer.
- The knee pad provides support for raising into standing and also to the user's legs during the transfer.
- Easily adjusted for access to wide chairs.
- The cross bar allows the patient to support themselves actively by pulling.

ReTurn 7500 SystemRoMedic (x1)



Summary: This is a multi-functional transportation platform that can facilitate patient transfer for those with mobility problems.

Age range: Adults

Max. Weight Limit: 150kg

Cleaning: wipe clean

Special features:

- Slip-resistant footplate
- Wheels lock for safety and security
- Adjustable leg support

Molift Raiser (x1)



Summary: The Molift Raiser Transfer System is a sit-to-stand patient-transfer system that has been designed to provide a safe and effective method of helping the user to their feet.

Cleaning: wipe clean

Max. Weight Limit: 170kg

Special features:

- Slip-resistant footplate
- Wheels lock for safety and security
- Adjustable leg support

Accessories:

- Safety strap with sliding sleeve, to assist user to get up into standing position.
- The Raiser Strap+ can be used by one or two carers when the user needs additional help to stand.

Ambiturn (x1)



Summary: The Ambiturn Sit-to-Stand Transfer Aid assists with pivoting transfers between different seat surfaces and can transport users short distances.

Special features: Footplate, handle, breaks, lower leg support.

Max. Weight Limit: 150 kg

Cleaning: Wipe clean

Mobile Hoist (x1)



Summary: Midi electric hoist features an emergency lift-and-lower function and a manual leg spreader, which allows the legs to be opened or closed to gain easy access around objects like beds, chairs, and wheelchairs.

Max. Weight Limit: 180kg (27 stone)

Cleaning: Wipe clean

Special Features:

Smart monitor feature which provides data to user and service engineer. Manual emergency descent function. Foot pedal leg opening mechanism with latch/lock functionality. On castors for easy manoeuvring. Detachable battery pack and desktop charger.

Room Use: Lillymere

Swivel Cushion (x2)




Summary: This swivel seat cushion assists with turning between bed wheelchair etc...

Max. Weight Limit: 150kg

Cleaning: remove cushion covers and wash in the washing machine

Comfort and Support

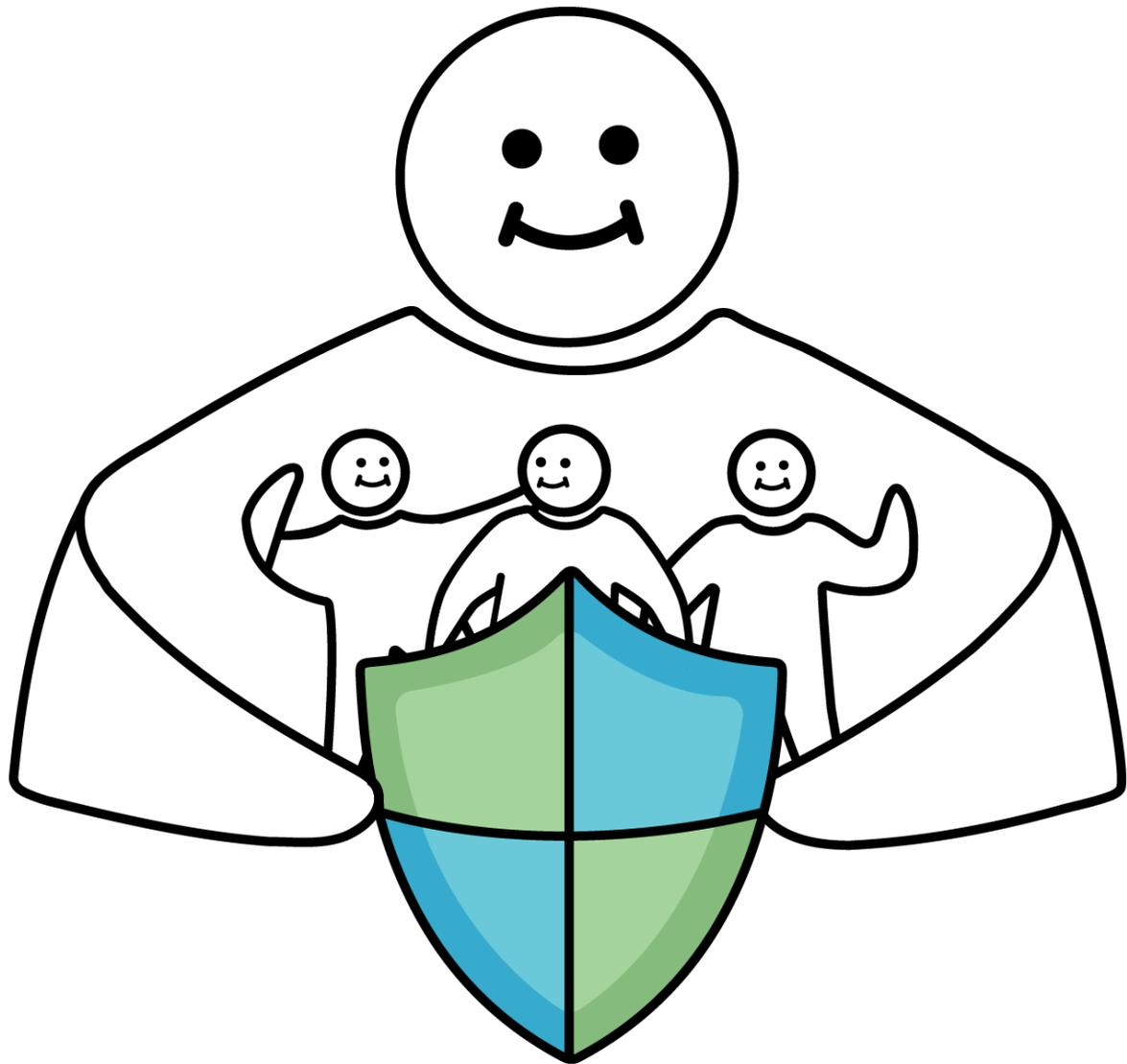

Soft-touch Tomato sitters (1x size 2 and 1x size 3)



Summary: Adaptive seating system that provides support in a variety of sitting positions

Age range: 1 – 10 years old (person dependant)

Important Info: Please do not lose the straps.

Cleaning: Wipe clean

Special features: Can be used in the pushchair, roughrider, canoe, dining chairs, and bikes.

Firefly GoTo (1x size 1)



Summary: Adaptive seating system that provides support in a variety of sitting positions

Age range: 1-3 years

Max. Weight Limit: 30kg

Important Info: Can be used with or without the base.

Cleaning: Washing Machine fabrics 30 deg C. Wipe down frame.

Special features:

Floor sitter has 26 recline positions, from upright (106°) to right back (140°). Height-adjustable headrest and laterals allow you to reposition the GoTo to fit child's size.

'P'Pods (Sizes 3-5.5)



Summary: The unique design of the 'P'POD provides the level of support required to stabilise posture, enabling clients with complex postural problems to sit comfortably.

Age range: 18 months -15 years, depending on size. We also have a custom large adult size one.

Important Info: Can be used for leisure in the lounge or bedroom as an alternative to wheelchair seating. Can also be used in canoes.

Cleaning: Wipe down

Special features: Can be used with or without head rest.

Accessories: Headrest and bolster beanbag.

Stabilo Cushions



Stabilo Grande, Stabilo Comforta, Stabilo comfortable plus duo

Summary:

Stabilo Grande is a mattress-sized cushion. Its versatility means that it may be used for lying, sitting or physiotherapy and rehabilitation. It is possible to adapt the same equipment for several users. The mattress is supplied with an electric Vacuum pump to quickly change and set the shape.

Stabilo comforta: The cushion can be used on its own to stabilise the back and hips, or together with the Stabilo Base cushion to create an anatomical seat.

Stabilo comfortable plus duo is a complex posture cushion which combines a corset back with a seat. Once shaped, it can be moved between a wheelchair, armchair or sofa.

Important Info: These cannot be booked out, but are available for use on activities.

Cleaning: Handwash , hose pipe in the bath or washing machine depending on size of cushion.

Special Features:

- Materials used are certified as antibacterial.
- Cushions are certified and non-combustible.
- Cushions are latex free

Electric Fans (x8)

Summary:

Freestanding, oscillating function.

Three speed setting and 90-degree oscillation.



Highchairs (x3)



Summary:

Has soft padded seat and a footrest for comfort. Folds flat.

Important Info:

Suitable from 6 months

Removable wipe-clean padded seat; 5-point harness with crotch post; Footrest;

Folds for easy storage;

Dimensions: D72 x W72 x H101cm

Cleaning: Wipe clean

Audio monitor (x3)

Video monitor (x4)



Summary: Video: 2.8" colour display with digital zoom. Infrared night vision allows you to see clearly with low light levels.

Important Info:

- Room temperature display.
- Notifies you of baby activity with audio alert.
- Up to 300m range and out of range indicator.
- Parent unit mains or battery powered.

Cleaning: Wipe clean

Mealtime support equipment



A selection of mealtime supports can be loaned. This must be returned to the kitchen. Dishwasher safe.

Deafguard fire alarm (x1)



Summary: Deafguard is a stand-alone, fire alarm listening device for deaf and hard of hearing people. Wire-free, battery-powered, acoustically-triggered bedside unit with high intensity strobe and vibrating pillow pad.

Cleaning: Wipe clean

Important Info:

The connected vibrator pad is placed under the pillow. When the Deafguard unit hears the fire alarm, the under-pillow pad vibrates, waking the sleeper, whilst the activated flashing strobe and LCD provides a complementary visual indicator in case the guest is in the room but not in bed. The result is an immediate and effective motion and visual trigger which signals that the fire alarm is sounding, encouraging the individual to leave the building and move to a place of safety.

Vibrating alarm clock facility.

Night lights (x8)



Summary:

Plug-in LED night light.

Cleaning: Wipe clean.

Lockable medicine fridge (x2)



Summary:

Having a secure lock on the medication fridge helps keep medications safe from unauthorised people.

This is a shared fridge and cannot be booked separately.

Location: 1 x Lodge building 1 x Acorn Building.

Cleaning: Wipe clean.

Smooth Talker



Summary: The Smooth Talker has 8 minutes of recording time, 5 levels and superior sound quality. It has a built-in rechargeable battery which will give many hours of use between charges and can be used while recharging.

Cleaning: Wipe clean

Special features:

Eight minutes of recording and playback time and five levels with adaptive recording time on each.

Built in rechargeable lithium battery - charge from a USB socket (can be used while charging).

Digital amplifier with twin speakers for superior sound quality.

Sequential, Random, Random (no repeats), Choice, Choice Progression, Converse, Auditory Prompt and Question/Answer playback modes.

Input socket for attaching a wired switch, or alternatively use a wireless switch such as the it-Switch.

Built-in symbol holder.

Large (125mm), low profile, switch activation area.

Stairgate/doorway safety barrier (x1)



Summary: FlexiGuard stairgates are a lightweight and portable alternative to conventional safety gates. 4 point pressure-fit safety barrier with quick release pressure-fit handle simultaneously locks all 4 pressure points into place.

Adjusts to fit openings between 71cm and 92cm.

Room use: Fits all rooms in the Lodge. Does not fit doors in Acorn House.

Cleaning: Wipe clean.

Door Alarms (X10 sets)



Summary: A door alarm is an alarm that is specifically set to go off when a door is opened. Door alarms are simple to install and are useful for carers of individuals who may leave their room. Door-open alarms activate when a door is opened by an individual and held open for a certain period of time. The door alarm warns the carer by emitting a chime, a loud siren sound or a flashing light.

Important Info: The alarm does not prevent individuals from wandering or leaving the room, but it does notify carers/parents if they do.

Cleaning: Wipe clean.

